



Registry Integration: STS Import Tool



STS Import Tool

- To use the STS import tool:
 - First use your STS vendor application to extract a harvest file containing MRNs (or SSNs if MRNs are not available).
 - The instructions for different vendors in subsequent slides.
 - Once you have a harvest file, use the MPOG Application Suite to input this data into your local MPOG database.
 - The data can then be PHI scrubbed and transferred to the central repository.

STS Harvest Schedule

Adult Cardiac Surgery Database - 2023

	Harvest Submission Window Close	Opt-Out Date	Includes Procedures Performed Through:	Report Posting	Comments
Harvest 1	2/10/2023	02/14/2023	12/31/2022	Spring 2023	Star Rating
Harvest 2	5/19/2023	05/23/2023	3/31/2023	Summer 2023	
Harvest 3	8/18/2023	8/22/2023	6/30/2023	Fall 2023	Star Rating
Harvest 4	11/10/2023	11/14/2023	9/30/2023	Winter 2023	

General Thoracic Surgery Database - 2023

	Harvest Submission Window Close	Opt-Out Date	Includes Procedures Performed Through:	Report Posting	Comments
Spring 2023	3/3/2023	3/7/2023	12/31/2022	Summer 2023	Star Rating
Fall 2023	9/8/2023	9/12/2023	6/30/2023	Winter 2023	Star Rating

Congenital Heart Surgery Database - 2023


	Harvest Submission Window Close	Opt-Out Date	Includes Procedures Performed Through:	Report Posting	Comments
Spring 2023	3/10/2023	3/14/2023	12/31/2022	Summer 2023	
Fall 2023	9/15/2023	9/19/2023	6/30/2023	Winter 2023	

STS Harvest Schedule


- Latest STS Harvest schedule available [here](#)

Vendor Application Instructions: ARMUS


1. File → Export → Select All
2. Pick Date Range: Chose field within Diagnoses and Procedures → Date of Surgery → add date range as filter
3. Fieldname Format: Short Name
4. Choice Format: Choose Harvest (codified values)
5. Field Export Order: Automatic
6. Export as a .dat file

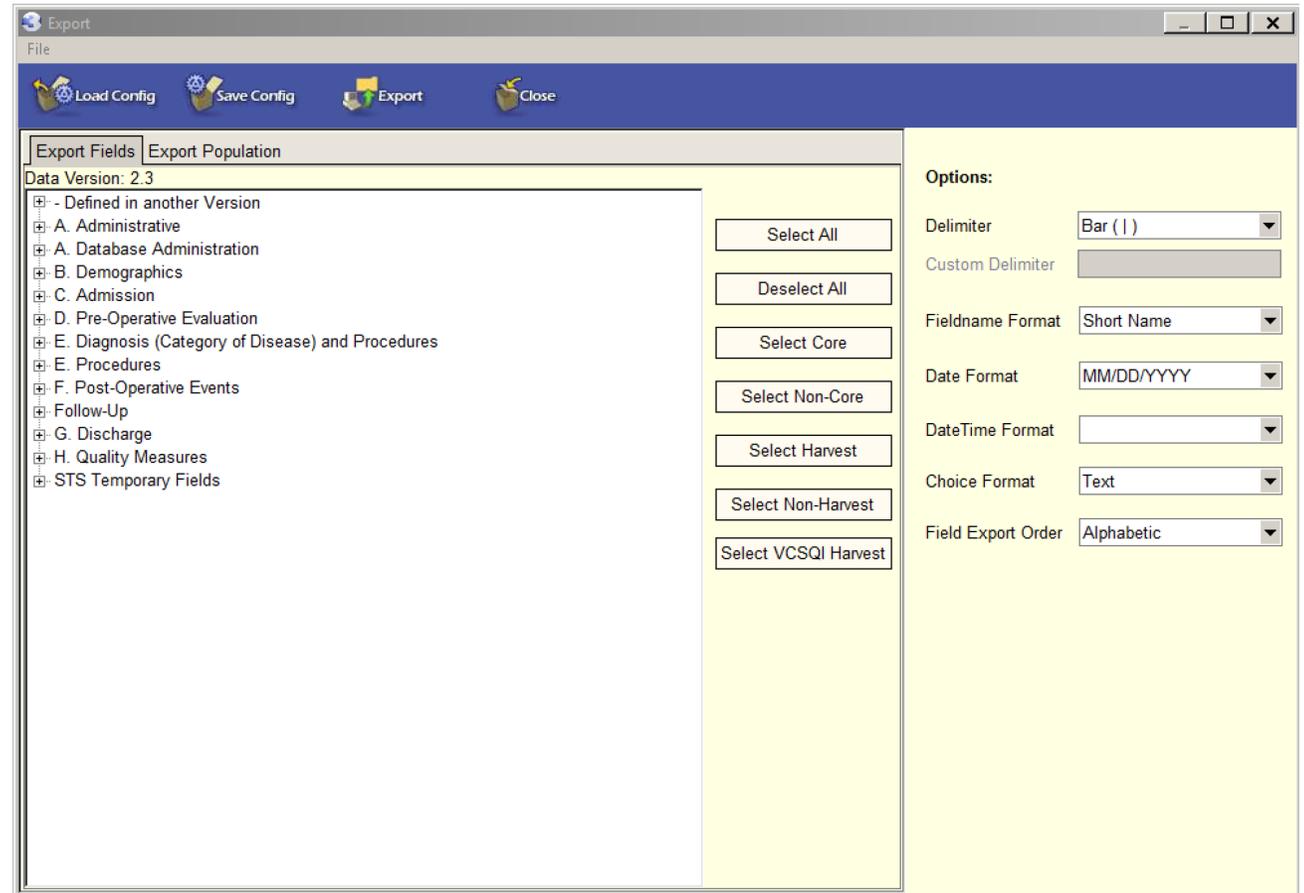

Vendor Application Instructions: ARMUS

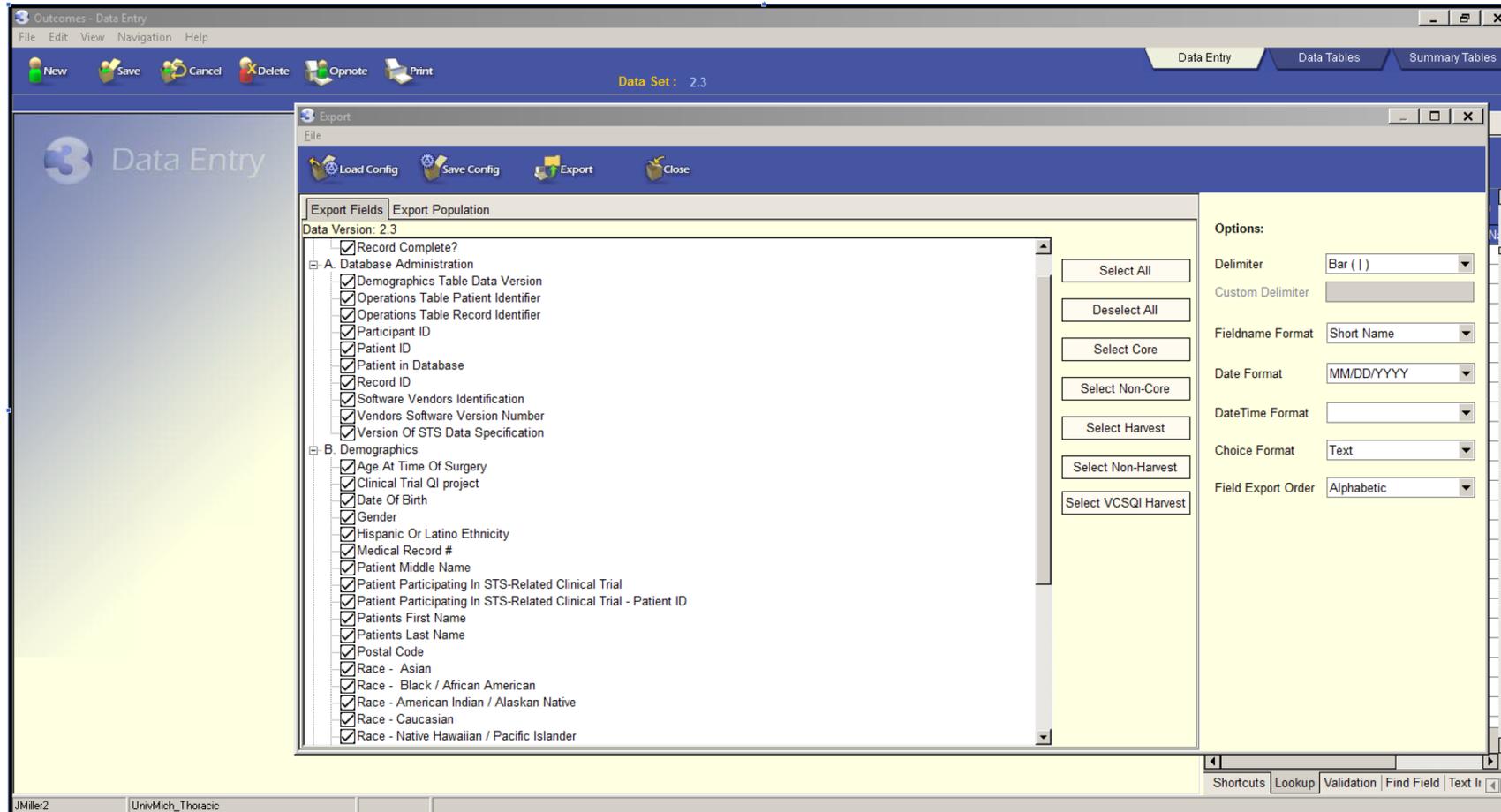


Vendor Application Instructions: ARMUS

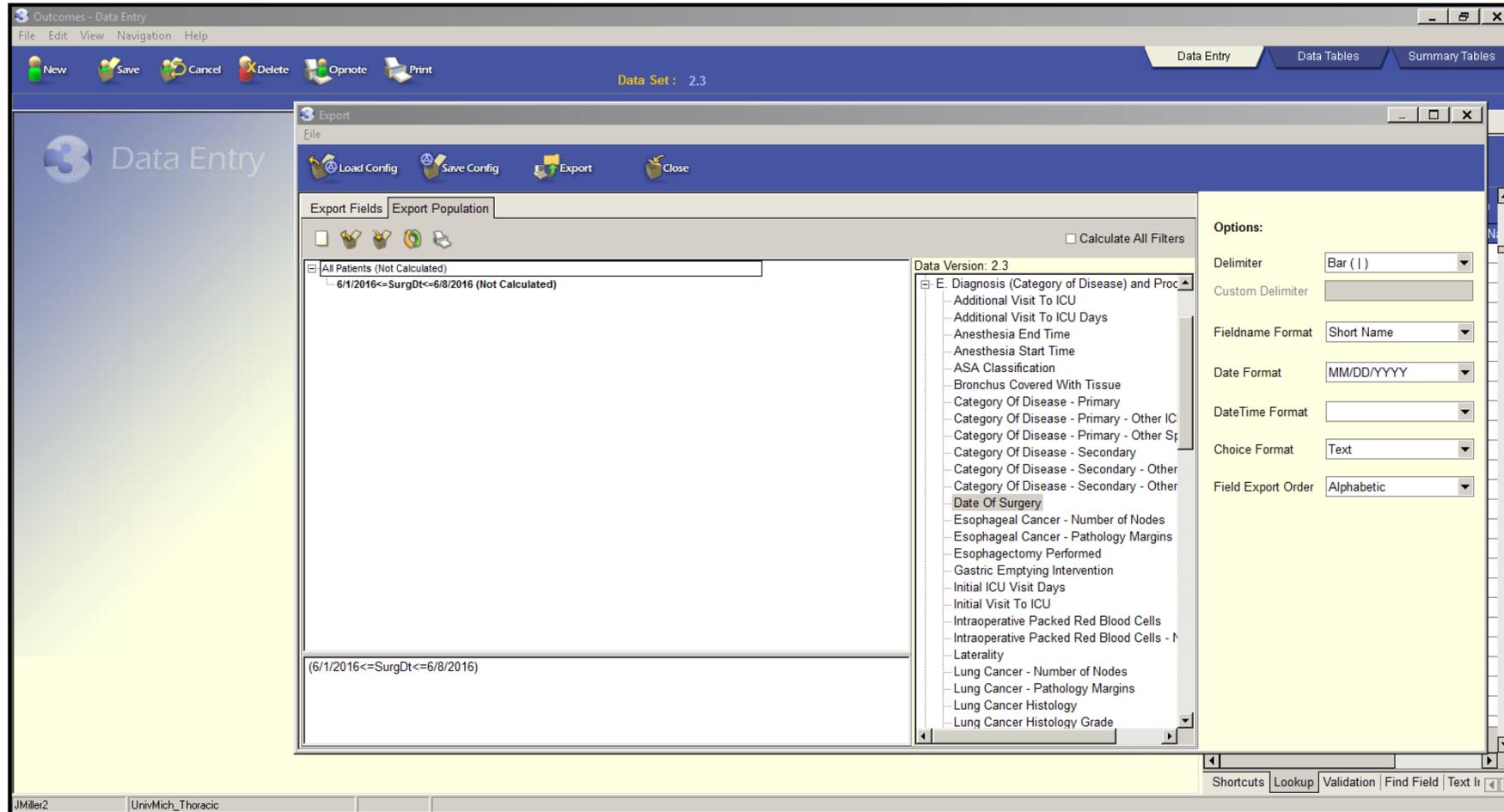


Vendor Application Instructions: LumedX

1. Select from the drop-down which STS version you would like to pull data for, then click the “play” button

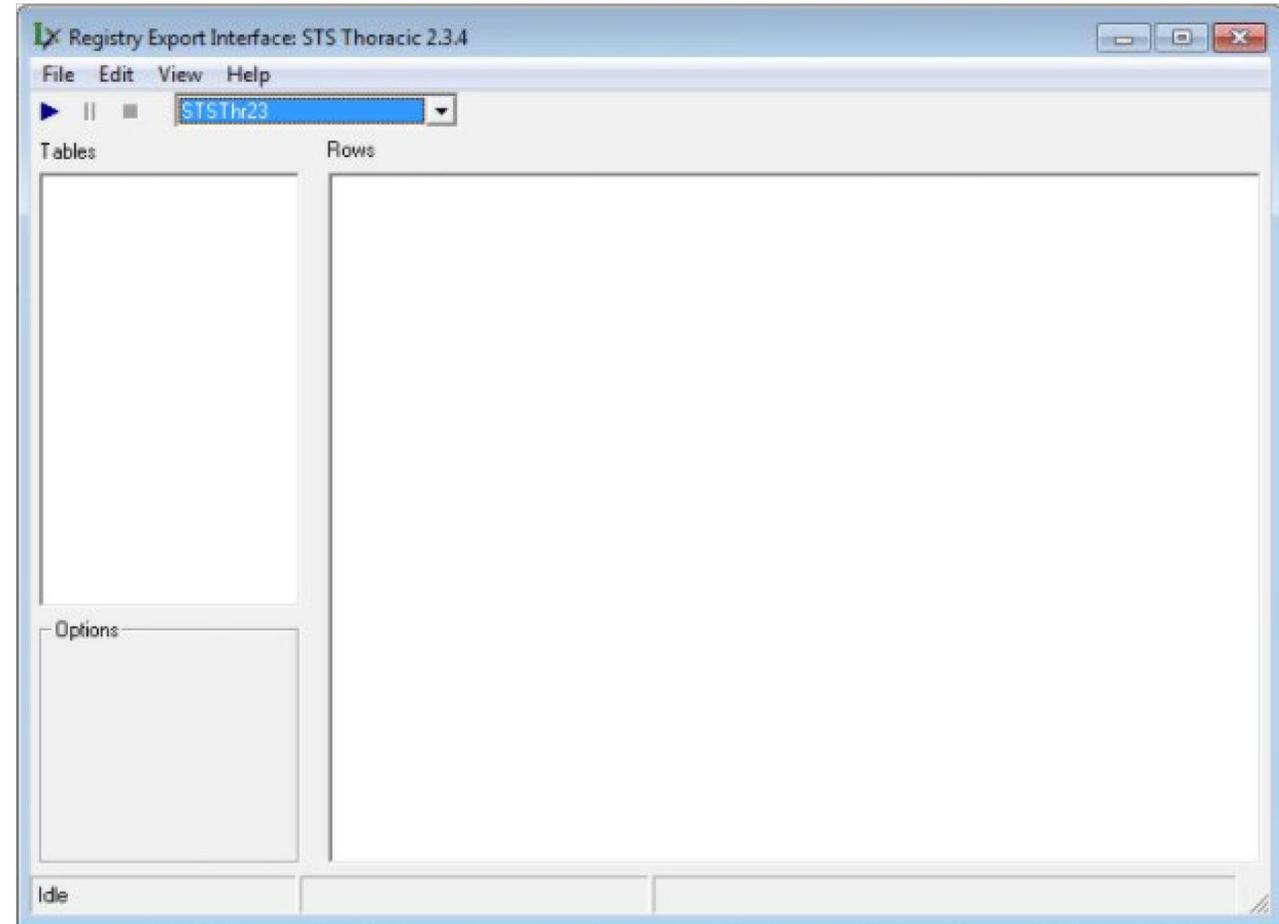


Vendor Application Instructions: LumedX

2. Enter your 5-digit participant ID into the box that pops up.

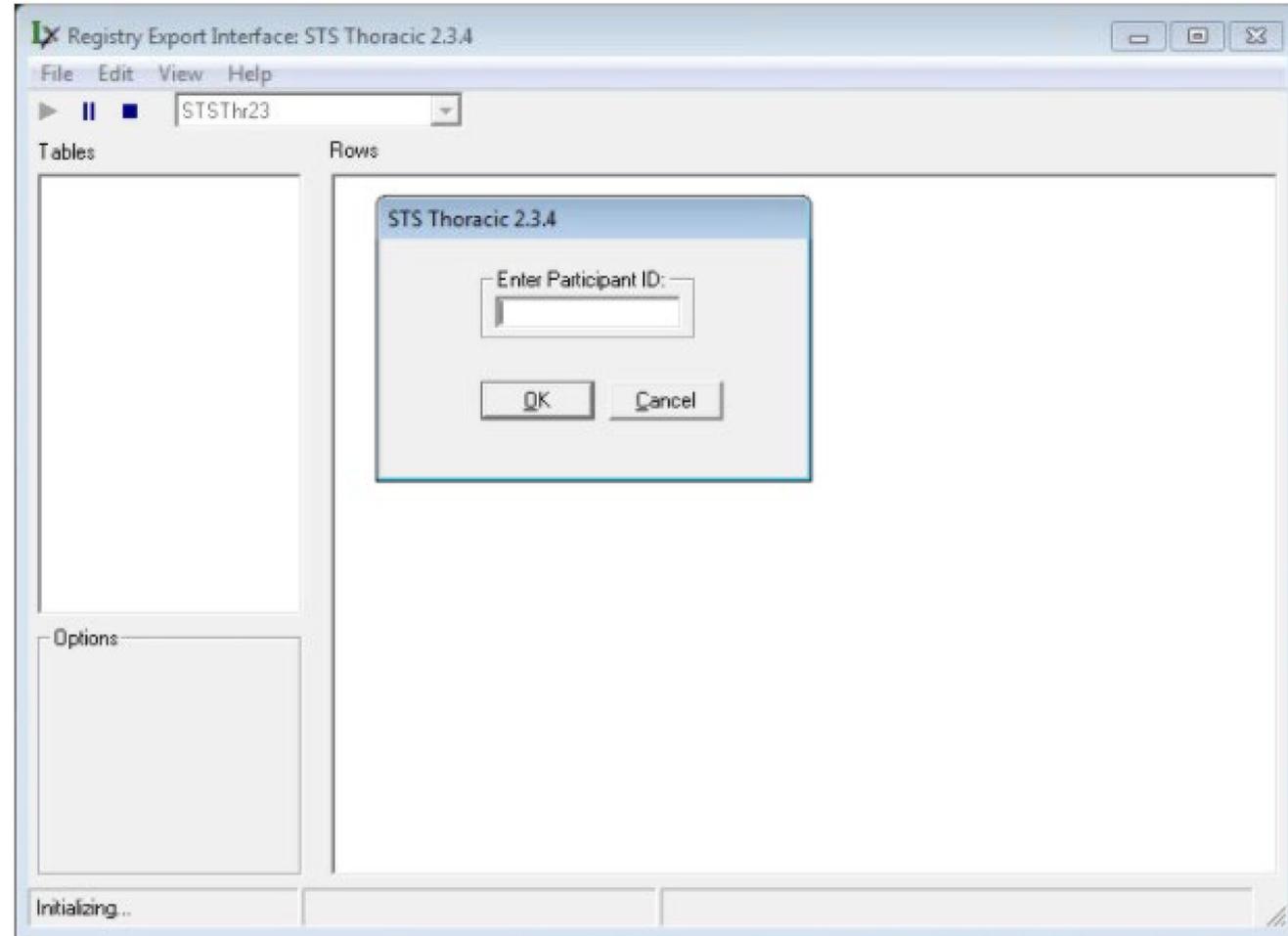


Vendor Application Instructions: LumedX

3. Specify the date range that you would like to extract data for, then click OK

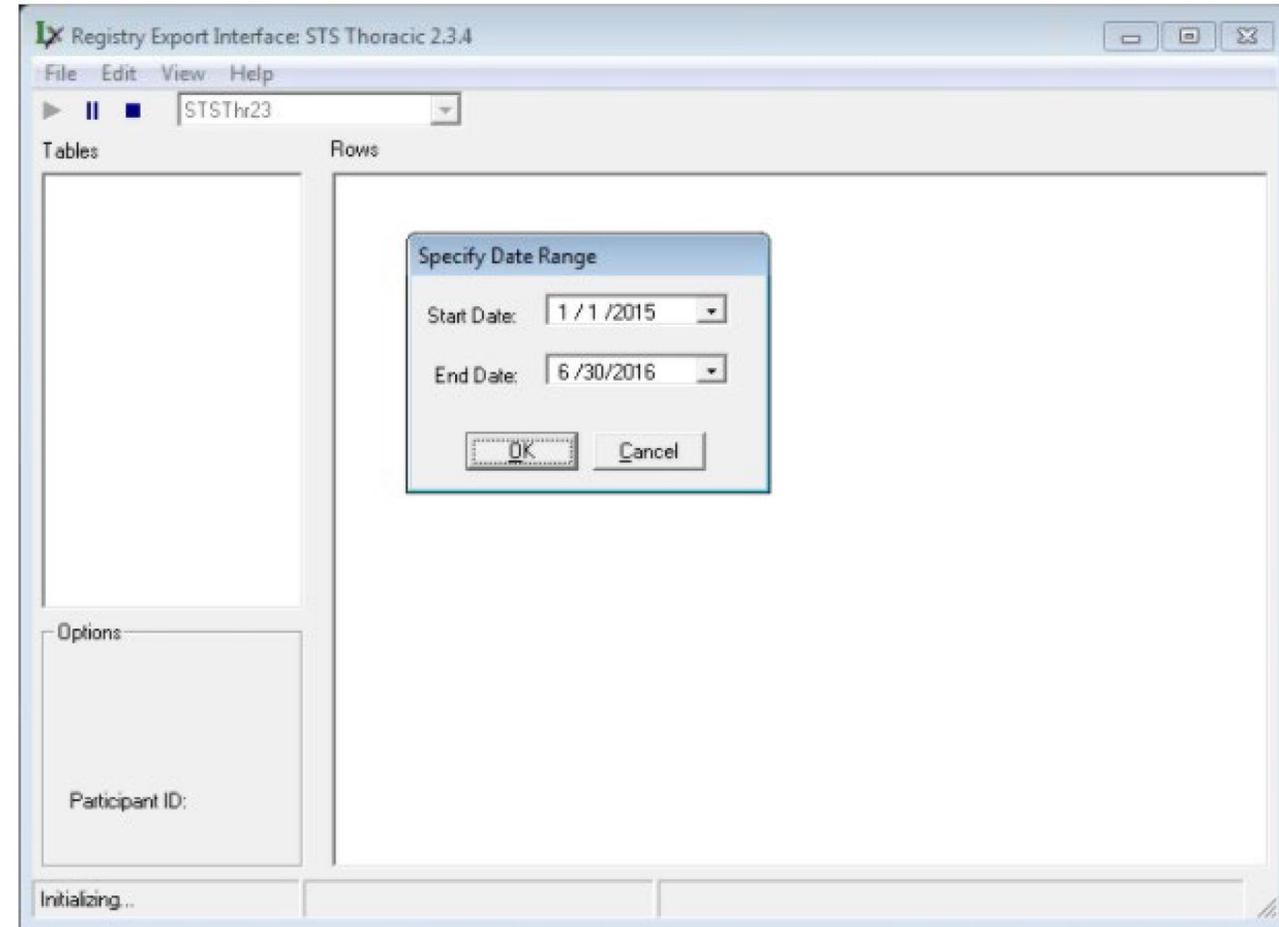


MPOG Import

1. Once you have your harvest file ready, open the MPOG Application Suite.
2. Open the STS Import application.

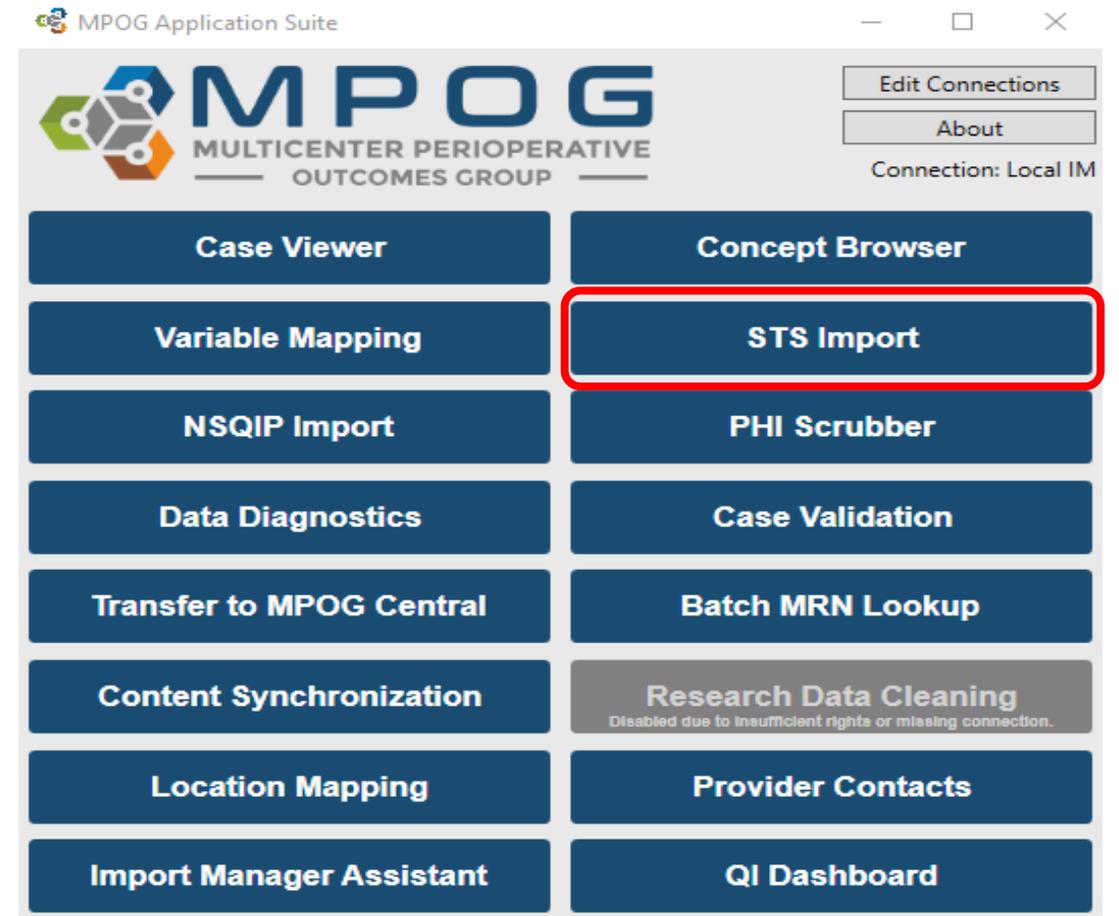

MPOG Import

3. The first screen will show you the import status (which months of data have been imported and what percentage of those STS records have been patient-matched locally). Click on the import tab to start an import of a harvest file

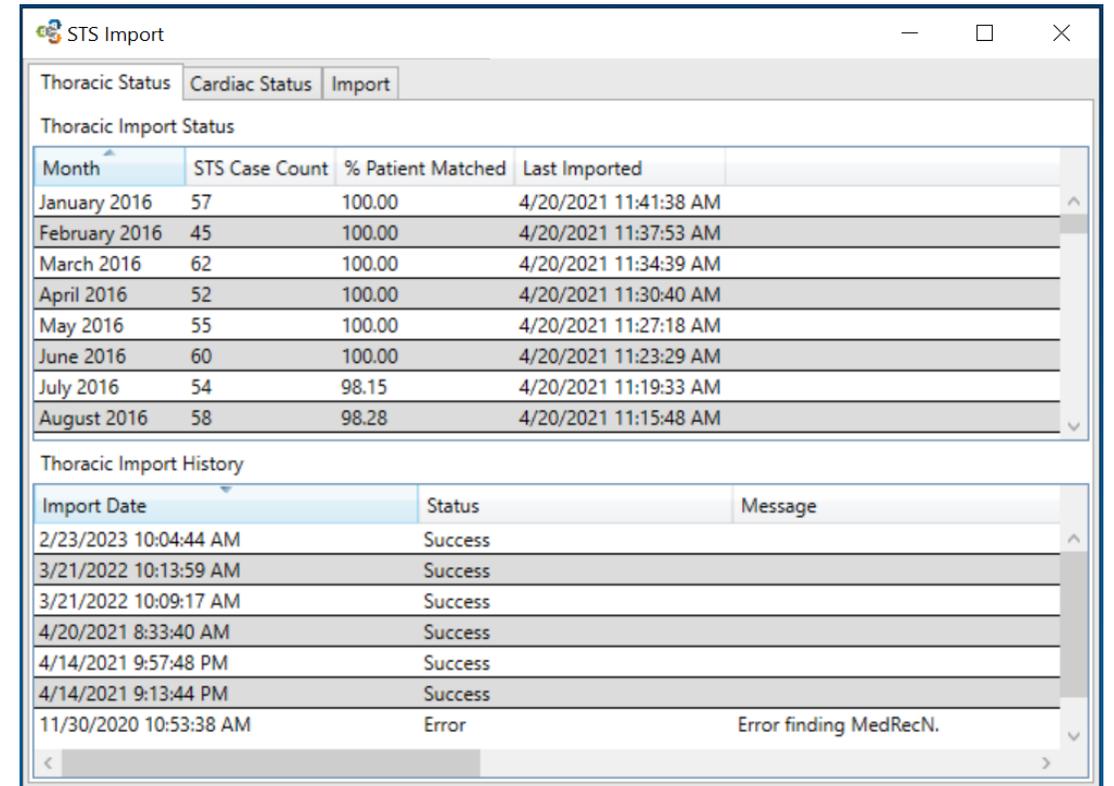

The screenshot shows the 'STS Import' application window. It has three tabs: 'Thoracic Status', 'Cardiac Status', and 'Import'. The 'Thoracic Status' tab is active, displaying two tables.

Thoracic Import Status


Month	STS Case Count	% Patient Matched	Last Imported
January 2016	57	100.00	4/20/2021 11:41:38 AM
February 2016	45	100.00	4/20/2021 11:37:53 AM
March 2016	62	100.00	4/20/2021 11:34:39 AM
April 2016	52	100.00	4/20/2021 11:30:40 AM
May 2016	55	100.00	4/20/2021 11:27:18 AM
June 2016	60	100.00	4/20/2021 11:23:29 AM
July 2016	54	98.15	4/20/2021 11:19:33 AM
August 2016	58	98.28	4/20/2021 11:15:48 AM

Thoracic Import History

Import Date	Status	Message
2/23/2023 10:04:44 AM	Success	
3/21/2022 10:13:59 AM	Success	
3/21/2022 10:09:17 AM	Success	
4/20/2021 8:33:40 AM	Success	
4/14/2021 9:57:48 PM	Success	
4/14/2021 9:13:44 PM	Success	
11/30/2020 10:53:38 AM	Error	Error finding MedRecN.

MPOG Import

4. Locate your STS harvest file using the “browse” button, then select if you are going to match on MRN or SSN (MRN preferred). Use the external matching CSV file option only in MRNs or SSNs are not populated directly within the harvest file

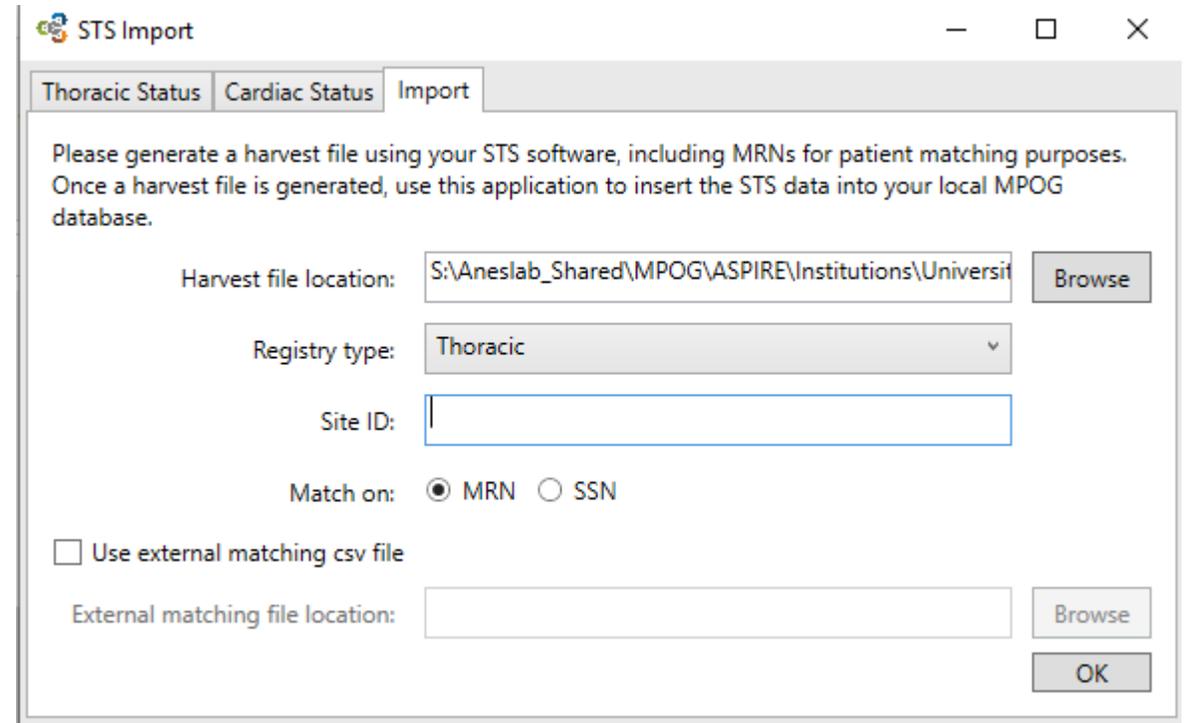


The screenshot shows the 'STS Import' application window with the 'Import' tab selected. The window contains the following fields and options:

- Harvest file location:** A text box containing 'S:\Aneslab_Shared\MPOG\ASPIRE\Institutions\Universit' and a 'Browse' button.
- Registry type:** A dropdown menu set to 'Thoracic'.
- Site ID:** An empty text box.
- Match on:** Radio buttons for 'MRN' (selected) and 'SSN'.
- Use external matching csv file**
- External matching file location:** An empty text box and a 'Browse' button.
- OK** button.

Instructions in the window: "Please generate a harvest file using your STS software, including MRNs for patient matching purposes. Once a harvest file is generated, use this application to insert the STS data into your local MPOG database."


MPOG Import

5. Click OK to start the import process. The process will import the STS data into the local MPOG database, perform patient matching, and set any PHI scrubber flags and transfer flags appropriately. After the process is complete, please run the PHI scrubber on “cases waiting for PHI scrub”

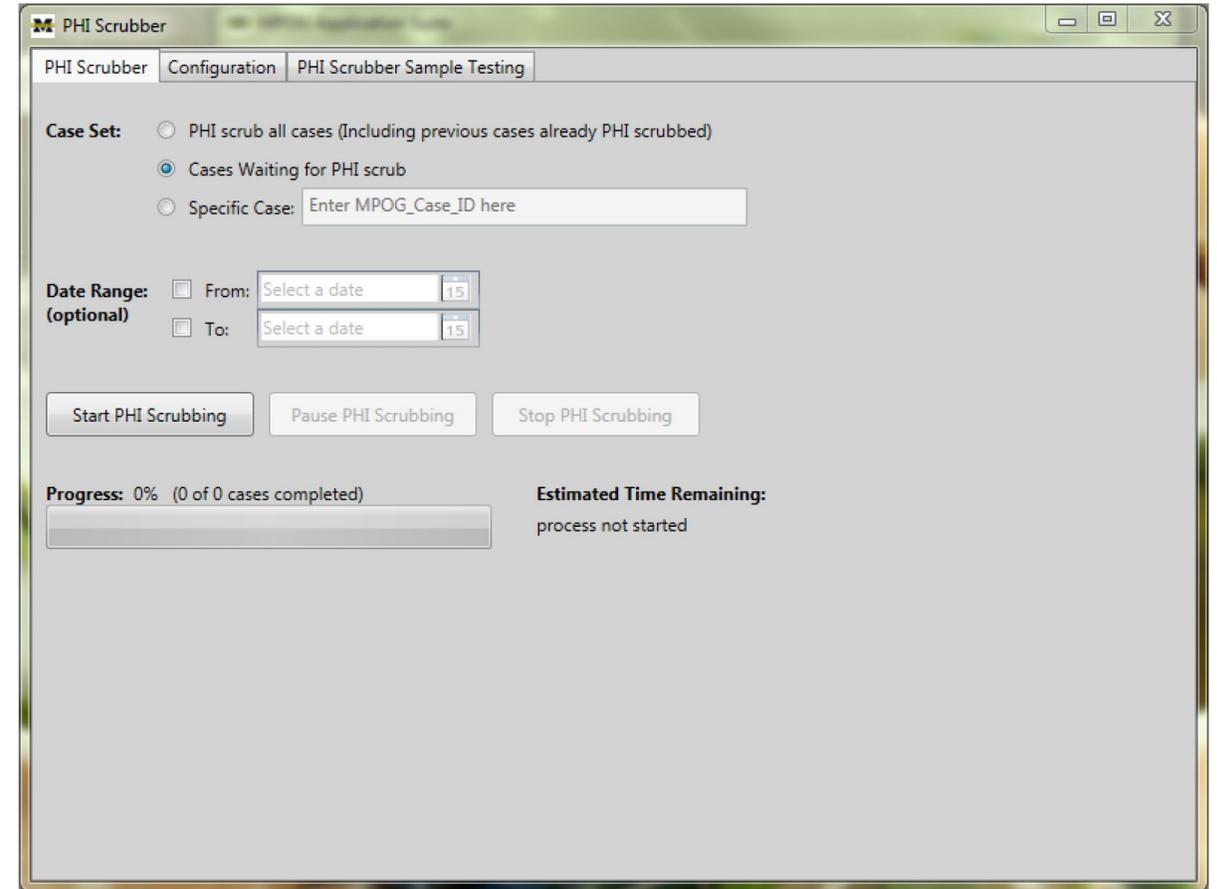

MPOG Import

6. To transfer the scrubbed STS records, please use the transfer application. Select the 'Historical' preset. You will also have to enter the dates that you want to transfer in the 'Specify Date Range' section.

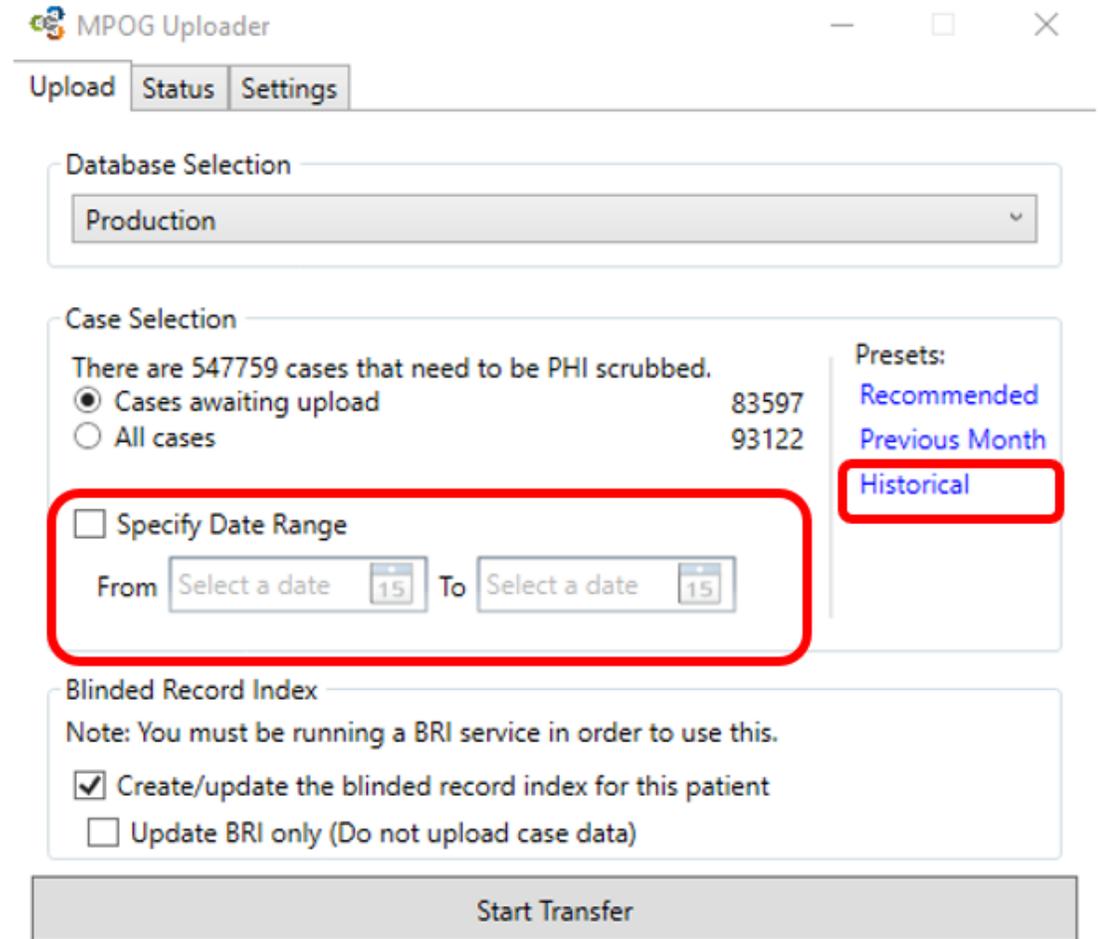

The screenshot shows the MPOG Uploader application window. It has three tabs: 'Upload', 'Status', and 'Settings'. The 'Upload' tab is active. The interface includes the following sections:

- Database Selection:** A dropdown menu showing 'Production'.
- Case Selection:** A section with the text 'There are 547759 cases that need to be PHI scrubbed.' It contains two radio buttons: 'Cases awaiting upload' (selected) with a count of 83597, and 'All cases' with a count of 93122. To the right, there is a 'Presets' section with three options: 'Recommended', 'Previous Month', and 'Historical'. The 'Historical' option is highlighted with a red box.
- Specify Date Range:** A section with a checkbox 'Specify Date Range' (unchecked). Below it, there are two date pickers: 'From' and 'To', both showing '15'.
- Blinded Record Index:** A section with the note 'Note: You must be running a BRI service in order to use this.' It contains two checkboxes: 'Create/update the blinded record index for this patient' (checked) and 'Update BRI only (Do not upload case data)' (unchecked).

At the bottom of the application is a large grey button labeled 'Start Transfer'.

Contact Info

- For questions, please contact support@mpog.zendesk.com for any technical assistance with the STS import application